**Waterpik® Ultra
Cordless Dental Water Jet**

**Hydrojet dentaire
Ultra Waterpik® sans fil**

**Irrigador dental
inalámbrico de
Waterpik® Ultra**

WELCOME

Congratulations! By selecting the **Waterpik® Ultra Cordless Dental Water Jet** you are now on your way to better gum health, and a healthier smile!

With a clinically proven combination of water pressure and pulsations, the **Waterpik® Dental Water Jet** removes harmful bacteria deep between teeth and below the gumline – where brushing alone won't reach.


We are so convinced that you will be happy with your new **Waterpik® Dental Water Jet**, that we back it up with a 14-day guarantee. If you aren't satisfied that your gums are healthier after using the **Waterpik® Dental Water Jet** for 14 consecutive days, feel free to return it to us for a full refund.

Thank you for choosing the **Waterpik® Ultra Cordless Dental Water Jet**.

Traditional Oral Care


Traditional brushing can't easily reach hidden food debris and harmful bacteria.

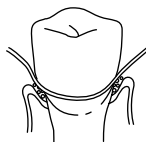Flossing has limited access below the gumline.

Oral Care with the **Waterpik® Dental Water Jet**

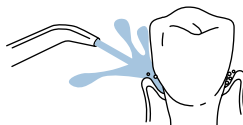The **Waterpik® Dental Water Jet** removes trapped food debris and harmful bacteria.

Important Safeguards	4
Product Description	6
Getting Started	8
Care of Your Waterpik® Ultra Cordless Dental Water Jet	13
Troubleshooting Guide	14
Limited Two-Year Warranty	15

IMPORTANT SAFEGUARDS

IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, the basic safety precautions below should always be followed.

READ ALL INSTRUCTIONS BEFORE USING.

DANGER:

To reduce the risk of electrocution:

- Do not handle charger with wet hands.
- Do not immerse in water or other liquid.
- Check the charger cord for damage before the first use and during the life of the product.
- Do not use while bathing.
- Do not place or store product where it can fall or be pulled into a tub, shower stall or sink.
- Do not reach for a product that has fallen into water. Unplug immediately.

WARNING:


To reduce the risk of burns, electric shock, fire, serious injury or damage to tissue:


- Do not use this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into any liquid. Contact Water Pik, Inc.
- Closely instruct and supervise children and individuals with special needs in the proper use of this product.
- Do not direct water under the tongue, into the ear, nose or other delicate areas. Without appropriate attachments recommended by Water Pik, Inc., this product is capable of producing pressures that may cause serious damage in these areas. See operating instructions for correct usage.
- Use this product only as indicated in these instructions or as recommended by your dental professional.
- Only use jet tips that are recommended by Water Pik, Inc.

- Fill reservoir with water or other dental professional recommended solutions only (i.e. mouthwash or prescribed solutions). A 1:1 ratio of mouthwash to water.
- Do not drop or insert any foreign object into any opening on the product.
- Keep charging unit and cord away from heated surfaces.
- Do not operate where oxygen or aerosol sprays are being used.
- Do not use mouthwashes/rinses that contain IODINE or Bleach.
- Do not use while wearing any oral jewelry. Remove prior to use.
- Do not use if you have an open wound on your tongue.
- If your physician advised you to receive antibiotic premedication before dental procedures, you should consult your dentist before using this instrument or any other oral hygiene aid.
- For household use only. Do not use outdoors.

SAVE THESE INSTRUCTIONS.

PRODUCT DESCRIPTION

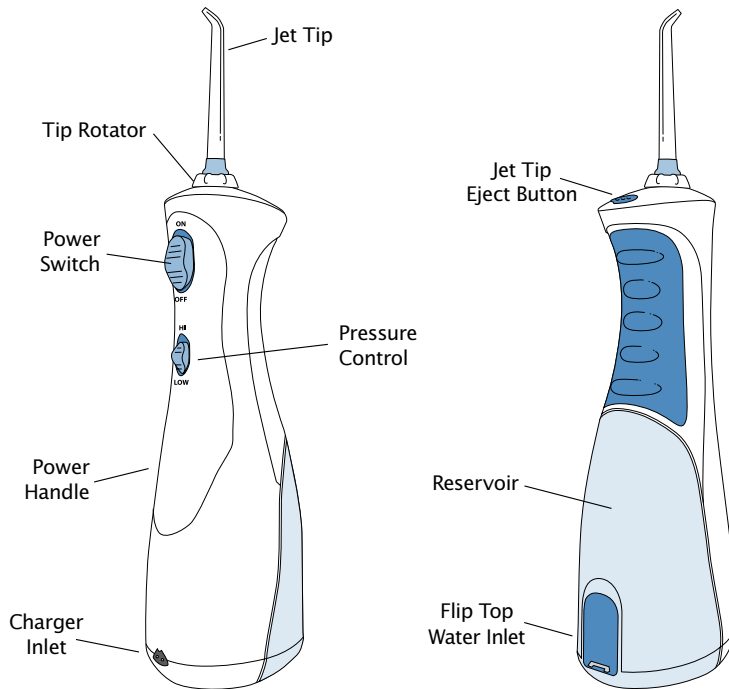
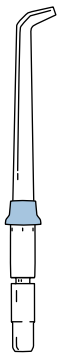
Standard
Jet Tip (1)


Tongue
Cleaner Tip (1)

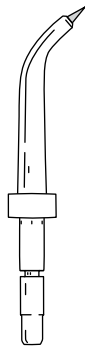Pik Pocket®
Tip (1)

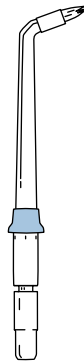Orthodontic
Tip (1)

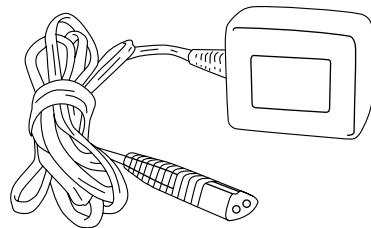
Charger

Replacement jet tips/attachments may be purchased directly from Water Pik, Inc.

Call 1-800-525-2774 (USA) or 1-888-226-3042 (Canada) or visit our web site at www.waterpik.com.


GETTING STARTED

Important! Turn Dental Water Jet off and charge for 16 hours before first use. Plug the cord into any standard electrical outlet. The Waterpik® Ultra Cordless Dental Water Jet will stay charged for approximately 1 week depending on use.

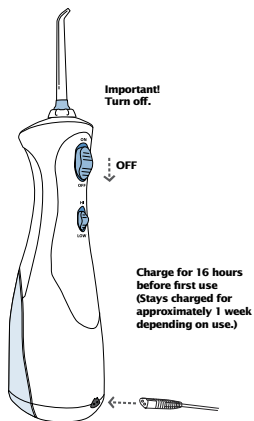
Filling the Reservoir

Lift flip top on reservoir and fill the reservoir with lukewarm water.

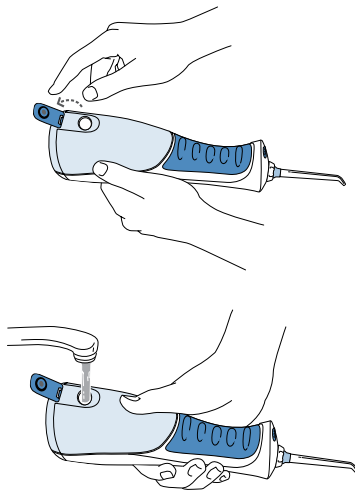Inserting and Removing Tips


Insert tip into the center of the knob at the top of the Dental Water Jet handle. The colored ring will be flush with the end of the knob if the tip is correctly locked into place.

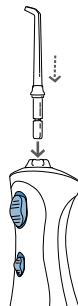
To remove the tip from the handle, press the jet tip eject button and pull tip from the handle.

Adjusting the Pressure Setting


You can adjust the pressure setting on the fly with your Ultra Cordless Dental Water Jet by simply sliding the pressure switch on the handle. Slide the pressure switch to “low” for gentle cleaning on sensitive gums. Slide the pressure switch to “high” for more powerful cleaning between teeth and below the gumline.

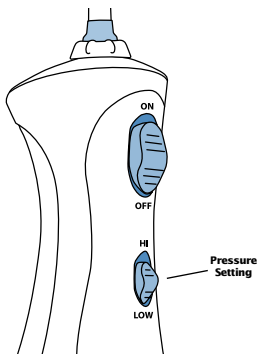Keep unit upright during use for best results.



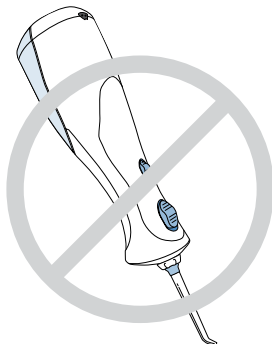
Direct the stream by turning the knob at the top of the handle.

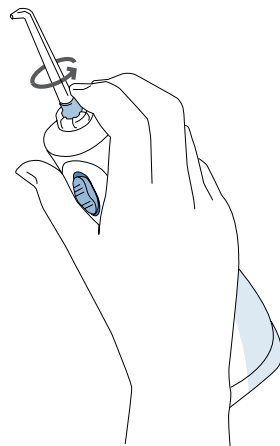
GETTING STARTED

Recommended Technique

Lean low over sink and place the jet tip in mouth. Aim the tip toward teeth, and turn the unit on.

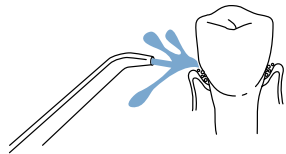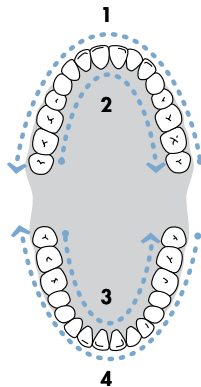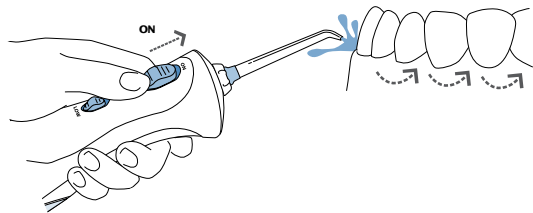
Direct the jet stream at a 90-degree angle to your gumline. **Slightly close lips to avoid splashing** but allow water to flow freely from mouth into the sink. For best results, start in the molar area (back teeth) working toward the front teeth. Follow the gumline and pause briefly between teeth. Continue until all areas around and between the teeth have been cleaned.

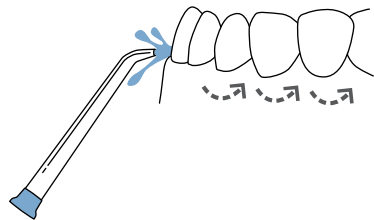
When Finished

Turn the unit off. Empty any liquid left in the reservoir.

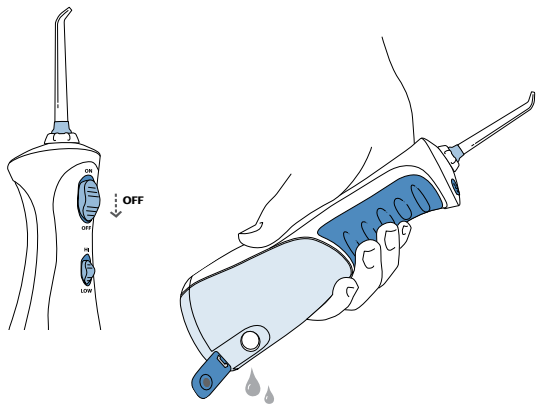If desired, the reservoir may be removed from the power handle by pulling the reservoir release latch while sliding the reservoir down towards the base of the unit.

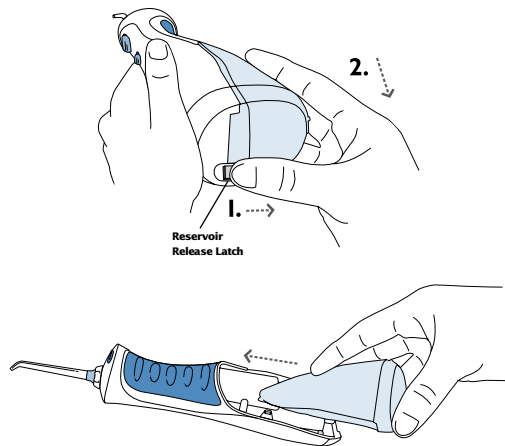To replace the reservoir on the power handle, simply slide it up towards the top of the unit. The reservoir will snap into place when fully engaged.



GETTING STARTED

SPECIAL TIPS

Pik Pocket® Tip

The Pik Pocket® Tip delivers water or anti-bacterial solutions deep into periodontal pockets.

To use the Pik Pocket® Tip, lean low over your sink and place the soft tip against a tooth at a 45-degree angle and gently place the tip under the gumline, into the pocket. Turn the unit on to the lowest pressure setting and clean by gently following gumline, inserting tip in between teeth. Continue tracing along the gumline until the reservoir is empty.

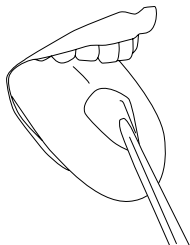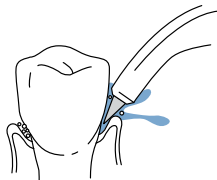
Tongue Cleaner Tip

Place the tongue cleaner on your tongue and start the Dental Water Jet – use the lowest pressure setting, increase as comfortable. Gently place the tongue cleaner in the center/midline of your tongue about half way back and with light pressure pull forward.

Orthodontic Tip

Start on the last molar (back tooth). Aim tip at a 90 degree angle to the tooth at the gumline. Gently glide tip along gumline, pausing briefly to lightly brush area between teeth and all around orthodontic bracket, before proceeding to the next tooth. Continue process until you have cleaned and irrigated the inside and outside of both the upper and lower teeth.

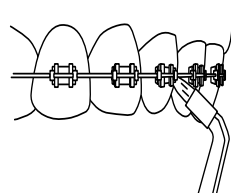For optimal results it is recommended that orthodontic tips are replaced after 6 months of use.

Replacement jet tips/attachments may be purchased directly from Water Pik, Inc.

Call 1-800-525-2774 (USA) or 1-888-226-3042 (Canada) or visit our web site at www.waterpik.com.

CARE OF YOUR WATERPIK® ULTRA CORDLESS DENTAL WATER JET

Never immerse the unit in water. Wipe the exterior when necessary with a soft cloth and mild non-abrasive cleanser. Reservoir may be washed with soapy water or in top rack of dishwasher.

Before cleaning, unplug from the electrical outlet.

Before exposing the Dental Water Jet to freezing temperatures, remove the reservoir and run the unit until completely empty.

Service Maintenance

Waterpik® Dental Water Jets have no consumer-serviceable electrical items and do not require routine service maintenance. Repair of this product should be performed only by an authorized service center or by Water Pik, Inc. For all your service needs go to www.waterpik.com and click on “product support” and “service locator” for your nearest service center. You may also call our toll-free customer service number at 1-800-525-2774 (USA) or 1-888-226-3042 (Canada).

Some replacement parts and accessories for the Ultra Cordless are also available for sale at the customer service numbers shown above.

Refer to the serial and model numbers in all correspondence. These numbers are located on the bottom of your unit.

TROUBLESHOOTING GUIDE

PROBLEM	CAUSE	CORRECTIVE ACTION
Jet Tips Leakage between the tip and the tip rotator.	1) Tip not fully engaged.	1) Remove tip and reinstall.
Dental System Reservoir leaks.	1) Liquid was left in reservoir. 2) Base not correctly engaged with reservoir. 3) Water inlet is open.	1) Always empty the reservoir after use. 2) Remove reservoir then reinstall as per instructions. 3) Close water inlet.
Inadequate pressure/No fluid flow.	1) Unit is upside down. 2) Insufficient liquid in reservoir. 3) Reservoir connection clogged. 4) Pressure setting on low.	1) Jet tip must be higher than reservoir for correct water flow. 2) Fill with liquid. 3) Fill reservoir with warm tap water and one cap of vinegar. Run until empty. Rinse with full reservoir of tap water. 4) Adjust pressure setting to high.
Unit doesn't start.	1) Battery not charged.	1) Be sure outlet is functioning; test with another appliance. Turn unit off. Plug in charger. Charge for 16 hours. 2) Be sure charger pins are clean and unobstructed.

Still have questions? Go to www.waterpik.com and click on "Dental Care" and the "Ask Ann" link.

If none of these corrective actions are successful, call the Water Pik, Inc. Customer Service Center. Toll-Free Number: 1-800-525-2774 (U.S.A.) 1-888-226-3042 (Canada).

LIMITED TWO-YEAR WARRANTY

Water Pik, Inc. warrants to the original consumer of this new product that it is free from defects in materials and workmanship for 2 years from the date of purchase. Consumer will be required to submit the original purchase receipt as proof of purchase date and if requested, the entire product, to support a warranty claim. Water Pik, Inc. will repair or replace at its discretion any part of the product, which in its opinion is defective, provided the product has not been abused, misused, altered or damaged after purchase, was used according to instructions, and was used only with accessories or consumable parts approved by Water Pik, Inc. Installation is the responsibility of the consumer and is not covered by the warranty. This limited warranty excludes accessories or consumable parts such as jet tips, etc.

This warranty gives you specific legal rights. You may also have other rights that vary, depending on the law in your place of residence.

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©2006 Water Pik, Inc.
Waterpik® and Waterpik® (stylized) and Pik Pocket® are registered trademarks of Water Pik, Inc.
treat yourself better™ is a trademark of Water Pik, Inc.
Waterpik®, Waterpik® (stylisé) et Pik Pocket® sont des marques déposées de Water Pik, Inc.
Waterpik®, Waterpik® (stylized) y Pik Pocket® son marcas comerciales registradas de Water Pik, Inc.
treat yourself better™ est une marque de commerce de Water Pik, Inc.